



70205 U.S. PTO  
08/803702  
02/21/97

Attorney Docket No.  
P-3639P1

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

ASSISTANT COMMISSIONER FOR PATENTS  
Washington, D.C. 20231

Sir: Transmitted herewith for filing is the patent application of:

Inventor(s): Vernon C. Maino and Maria Suni

For: Method for Detecting T Cell Response to Specific Antigens in Whole Blood

Also enclosed are:

- ☒ 3 sheet(s) of drawing(s). ☐ Formal ☒ Informal
- ☐ A Declaration and Power of Attorney ☐ Executed ☐ Unexecuted
- ☐ An assignment of the invention to Becton Dickinson and Company.
- ☐ A certified copy of a application.
- ☐ An Associate Power of Attorney to:
- ☐ An Information Disclosure Statement under 37 CFR §1.97.
- ☐ Copies of documents cited in the Information Disclosure Statement.
- ☐ A Sequence Listing under 37 CFR §1.821(c).
- ☐ A computer-readable form of the Sequence Listing under 37 CFR §1.821(e).
- ☐ A statement under 37 CFR §1.821(f).

CLAIMS AS FILED				
(1) FOR	(2) NUMBER FILED	(3) NUMBER EXTRA	(4) RATE	(5) BASIC FEE \$770.00
TOTAL CLAIMS	18 - 20	0	X \$22 =	\$0.00
INDEPENDENT CLAIMS	1 - 3	0	X \$80 =	\$0.00
MULTIPLE DEPENDENT CLAIMS	IF ONE OR MORE ADD \$250 TO COLUMN (5)		\$260	
	TOTAL FILING FEE →			\$770.00

- ☒ Please charge my Deposit Account No. 02-1666 in the amount of \$770.00. Triplicate copies of this sheet are enclosed.

"EXPRESS MAIL" LABEL NO. TB617117545US

DATE OF DEPOSIT February 21, 1997

I HEREBY CERTIFY THAT THIS PAPER OR FEE IS BEING DEPOSITED WITH THE UNITED STATES POSTAL SERVICE "EXPRESS MAIL POST OFFICE TO ADDRESSEE" SERVICE UNDER 37 CFR 1.10 ON THE DATE INDICATED ABOVE AND IS ADDRESSED TO ASSISTANT COMMISSIONER FOR PATENTS, WASHINGTON, D.C. 20231


Charles Miller  
NAME OF PERSON MAILING PAPER OR FEE

Charles Miller  
SIGNATURE OF PERSON MAILING PAPER OR FEE

☒ The Assistant Commissioner is hereby authorized to charge any additional fees which may be required, or credit any overpayment to Deposit Account No. 02-1666. Triplicate copies of this sheet are enclosed.

☐ A check in the amount of \$770.00 to cover the filing fee is enclosed.

Becton Dickinson and Company  
1 Becton Drive  
Franklin Lakes, NJ 07417-1880  
(201) 847-6722

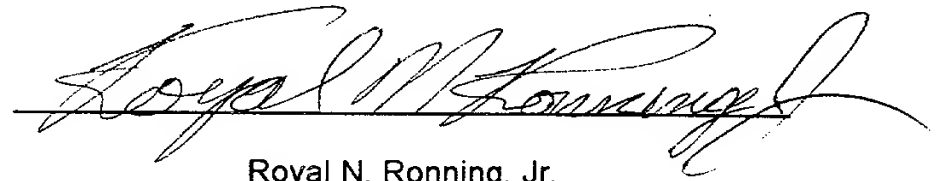

Royal N. Ronning, Jr.  
Attorney/Agent for Applicant(s)  
Registration No. 32,529

08/803702



02/21/97

Attorney Docket No.

P-3639P1



## IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

ASSISTANT COMMISSIONER FOR PATENTS  
Washington, D.C. 20231

Sir: Transmitted herewith for filing is the patent application of:

Inventor(s): Vernon C. Maino and Maria Suni

For: Method for Detecting T Cell Response to Specific Antigens in Whole Blood

Also enclosed are:

- ☒ 3 sheet(s) of drawing(s). ☐ Formal ☒ Informal
- ☐ A Declaration and Power of Attorney ☐ Executed ☐ Unexecuted
- ☐ An assignment of the invention to Becton Dickinson and Company.
- ☐ A certified copy of a application.
- ☐ An Associate Power of Attorney to:
- ☐ An Information Disclosure Statement under 37 CFR §1.97.
- ☐ Copies of documents cited in the Information Disclosure Statement.
- ☐ A Sequence Listing under 37 CFR §1.821(c).
- ☐ A computer-readable form of the Sequence Listing under 37 CFR §1.821(e).
- ☐ A statement under 37 CFR §1.821(f).

## CLAIMS AS FILED

(1) FOR	(2) NUMBER FILED	(3) NUMBER EXTRA	(4) RATE	(5) BASIC FEE \$770.00
TOTAL CLAIMS	18 - 20	0	X \$22 =	\$0.00
INDEPENDENT CLAIMS	1 - 3	0	X \$80 =	\$0.00
MULTIPLE DEPENDENT CLAIMS	IF ONE OR MORE ADD \$250 TO COLUMN (5)		\$260	
	TOTAL FILING FEE →			\$770.00

- ☒ Please charge my Deposit Account No. 02-1666 in the amount of \$770.00. Triplicate copies of this sheet are enclosed.

"EXPRESS MAIL" LABEL NO. TB617117545US

DATE OF DEPOSIT February 21, 1997

I HEREBY CERTIFY THAT THIS PAPER OR FEE IS BEING DEPOSITED WITH THE UNITED STATES POSTAL SERVICE "EXPRESS MAIL POST OFFICE TO ADDRESSEE" SERVICE UNDER 37 CFR 1.10 ON THE DATE INDICATED ABOVE AND IS ADDRESSED TO ASSISTANT COMMISSIONER FOR PATENTS, WASHINGTON, D.C. 20231

Charles Miller

NAME OF PERSON MAILING PAPER OR FEE


Charles Miller

SIGNATURE OF PERSON MAILING PAPER OR FEE


☒ The Assistant Commissioner is hereby authorized to charge any additional fees which may be required, or credit any overpayment to Deposit Account No. 02-1666. Triplicate copies of this sheet are enclosed.

☐ A check in the amount of \$770.00 to cover the filing fee is enclosed.

**Becton Dickinson and Company**  
1 Becton Drive  
Franklin Lakes, NJ 07417-1880  
(201) 847-6722

  
\_\_\_\_\_  
Royal N. Ropping, Jr.

Royal N. Ronning, Jr.  
Attorney/Agent for Applicant(s)  
Registration No. 32,529

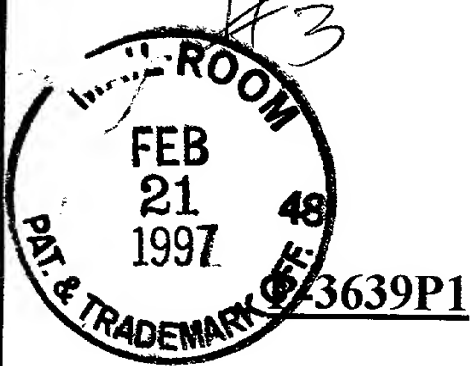

A

## METHOD FOR DETECTING T CELL RESPONSE TO SPECIFIC ANTIGENS IN WHOLE BLOOD

5

This application is a continuation-in-part of copending application Serial Number 08/760,447 filed December 6, 1996.

### BACKGROUND OF INVENTION

10

Despite widespread agreement that the basic unit of function of the T cell arm of the cellular immune system is the collection of memory/effector T cell clones responsive to a specific antigen(Ag), there is very little understanding of the development and homeostasis of Ag-specific human memory/effector T cells in either health or disease. In large part, this lack of definitive information is methodological; i.e. the simple lack of tools to reliably quantitate and functionally evaluate Ag-specific T cells in the human. For many years, assays of antigen specific T cell function relied on proliferation of cultured PBMC as the measure of the responding T cells.

20

New knowledge of immunological processes in the past decade, however, have highlighted serious problems with this approach. Proliferative responses require days to manifest, and it is now well documented that many, if not most, memory/effector T cells respond to Ag with an effector response (e.g. cytokine release or cytotoxicity) and then succumb to activation-induced apoptosis prior to any possible proliferative response. Those cells that survive to proliferate are thus a highly variable subset (depending on culture conditions) of the overall Ag-

25

responsive cohort. This issue is even more critical for the evaluation of T function in diseases such as AIDS, as it has been clearly documented that memory/effector T cells from HIV+ individuals display an increased susceptibility to apoptosis, and indeed, show enhanced Ag-induced apoptosis *in vitro*. Similarly T cell responses to antigen in other diseases may be characterized by varying ratios of cells undergoing proliferation, cytokine expression, apoptosis, or cytotoxic T cell differentiation.

Recently, other assays of T cell function have also included cytokine release endpoints, a defining function of memory/effector T cells. While this addition to the measurement of proliferative responses represents a clear advance, in terms of a more comprehensive measure of T cell activation, these assays do not define the functional cell types that participate in these responses. Further, most of these analyses have relied on ELISA-based quantitation of expressed levels of cytokine in culture supernatants. These procedures yield only bulk measurements; individual cell types participating in the response are not evaluated.

To assess which cell types are expressing selected cytokines, subset purification or single cell-based methodologies such as the ELISPOT and cytokine-release limiting dilution assays have been employed. However, these assays are extremely laborious -- especially if set up to measure the multiple, potentially relevant cytokines which may be produced by functionally heterogeneous Ag-specific memory/effector T cells -- and therefore are not amenable to routine use in the clinical laboratory. Moreover, the sensitivity of these assays may be limited by (1) conditions which do not mimic *in vivo* cellular environments which may be more clinically relevant, (2) T cell stimulation conditions that must necessarily compromise between conditions optimal for T cell activation and those compatible

with the single cell detection strategy of the assay (e.g. limiting dilution conditions or spot counting on nitrocellulose covered culture wells), and (3) vulnerability to some loss of Ag-reactive cells from activation-induced apoptosis (as these assays require from 1-2 to many days of *in vitro* culture, and activation-induced apoptosis, as measured by DNA strand breaks, can initiate in as few as 12-19 hours).

It has been previously shown that T cell cytokine responses occur in cells expressing the CD69 antigen. A TH1, or type I response, for example is typically identified as the percentage of CD3+ T cells co-expressing CD69 and  $\gamma$ -IFN, IL-2, or TNF- $\alpha$ . A TH2, type II, response can be characterized as the frequency of CD3+ T cells which co-express CD69, and IL-4 or IL-13. The distribution of individual cytokines of a type I or type II response may also be descriptive of a specific kind of cellular response to antigen. For example unique T cell subsets may express TNF- $\alpha$  and  $\gamma$ -IFN in response to different antigens and have been referred to as TH1 cells.

### **SUMMARY OF INVENTION**

This invention comprises a novel approach to the assessment of Ag-specific T cells that quantitates and characterizes these cells with unprecedented clarity, and importantly, because it is performed in whole blood, is amenable to routine use in the clinical immunology laboratory. In order to evaluate clinical samples for the types of antigen specific cellular responses that reflect *in vivo* conditions, we have developed an improved flow cytometric intracellular cytokine assay in whole blood

that can simultaneously measure multiple T cell subsets expressing multiple cytokines from a single whole blood culture. Evaluation of whole blood antigen specific cytokine responses has the important advantage of assessing T cell activation in the presence of ALL types of MHC autologous antigen presenting cells present in the native sample. It also has the advantage of enabling a culture system (whole blood) which can reflect effects of systemic environments (i.e. drug augmentation or suppression) on T cell responses to specific stimuli including antigen. by either culturing in the presence of such drug or analyzing the blood of a human or animal receiving such drug.

At its simplest, the methodology involves a step process, which involves culturing with the antigen specific stimulus and analyzing an aliquot of the cultured sample for expression of one or more intracellular cytokines and/or early activation antigens in combination with one or more T-cell markers, optionally with the lysing of the red blood cells and washing to remove debris.

The method described herein is composed of two parts. The first process relates to the cell culture process for optimal antigen presentation in MHC directed T cell responses (primarily but not restricted to memory CD4+ T cell responses). The second part relates to the analytical (flow cytometry) procedure which enables detection of antigen specific responding T cells based on the expression of an external activation antigen and an intracellular functional cytokine.

More specifically, this invention provides an assay protocol using whole blood for the rapid (generally less than 24 hours, preferably less than 6 hours), highly efficient, Ag-specific activation of secretion-inhibited CD4+ (memory/effector) T cells, followed by the quantitation and characterization of these



Ag-specific T cells using multiparameter flow cytometric, immuno-fluorescent detection of one or more intracellular cytokines (including IL-2, IL-4,  $\gamma$ -IFN and TNF-  $\alpha$  ) and an early activation antigen (such as CD69), in combination with one or more T cell subset-defining phenotypic markers (such as CD3 or CD4). It is further disclosed that the antigen specific response can be more easily detected when the antigen stimulation is provided in conjunction with costimulation of surface antigens involved with accessory cell surface molecules, such as CD28, CD40, VLA-4, and other such specificities known in the art. Costimuli can be either antibodies or ligands binding to these antigens. A preferred costimulus is CD28. It is also determined that T cell responses measured using the protocol defined in this disclosure can be shown to be sensitive to drugs which have been demonstrated to augment or suppress cellular responses (e.g. exogenous cytokines, cyclosporine-A, herbimycin-A).

Further, it has been found intracellular cytokine detection is enhanced when an agent which blocks the secretion of such intracellular cytokines is added during the activation (4h) period of incubation. A preferred agent is Brefeldin A.

T cells are activated by antigen by autologous antigen presenting cells (APC) (B cells, monocytes, or dendritic cells) which process protein antigen and present resultant peptide determinants on cell surface MHC molecules. We have determined that T cell/APC interactions result in inflammatory reactions that promote enhanced adhesion of APC/T cell aggregates. To detect cognate T-cell responses to specific antigen in suspension for flow cytometric analysis, cation chelators are employed to release reactive T cells into the analysis buffer. EGTA or EDTA are preferred compounds for this purpose.

A preferred methodology of this invention is presented below:

### **Antigen activation and culture of whole blood**

1. Aliquot 1 ml of heparinized whole blood into tissue culture tubes
2. For enhanced frequency of antigen specific T cell response, add 4µg/ml of CD28 (1µg/10<sup>6</sup> cells), incubate 15 min. at RT
3. Add 60µl/ml of antigen (typically 1 to 10µg final per 1 ml culture), incubate for 1h at 37C, 7% CO<sub>2</sub> (slanted at 5° angle)
4. Add 10µg/ml of BFA, incubate an additional 5h, with or without agitation
5. Add 100 µl of 20mM Na<sub>2</sub>EDTA, final concentration 2mM, vortex, incubate for 15 min. @ RT, wash with PBS with vigorous pippetting.
6. Lyse and fix RBC; e.g. with FACSlysing solution; wash in PBS
7. Proceed to “Staining and detection” , or:
8. Freeze in 10% DMSO/ 1% BSA/ PBS at -70C

### **Staining and detection of antigen responsive T cells after culture in whole blood.**

1. Thaw cells rapidly in 37C water bath
2. Wash 2X with cold PBS (at 500 x g) for 5 min.

3. Aliquot into FACS tubes, wash with FACS wash buffer
4. Permeablize cells with 500 µl of FACS permeabilization solution™ for 10 min. at RT

5

5. Wash with FACS was buffer

6. Stain with MAbs defining intracellular or surface components (e.g. CD3, CD4, cytokines, Brefeldin blocked CD69). Typically CD69 is conjugated to phycoerythrin and anti-cytokine MAbs (IL-2, γ-interferon, and TNF-α) are conjugated with FITC (see examples).

7. Wash and fix with 1% PFA

8. Analyze on FACScan (fluorescence trigger on T cell subsets using CD4 PerCP, or CD3 PerCP collect a minimum of 50,000 CD4+ T cell events.

It is to be understood that the above procedures may be modified as conditions dictate, and that a variety of reagents may be used in the permeabling and washing steps as would be apparent to the artisan.

For primary cell surface staining of whole blood cultures the procedure is modified as follows: This modification may not be necessary if cell surface epitopes

can survive permeabilization procedure. Examples of cell surface-defining reagents that may also be used in intracellular protocols include Leu-4 (CD3), Leu-3a (CD4), and Leu-2a (CD8).

(1) Stain cell surface antigens on T cells using titrated concentrations of appropriate monoclonal antibody (e.g. CD45RO) conjugated with fluorochrome; wash and fix with FACSlyse™.

(2) Permeabilize cells exposed to lysing solution with FACS permeabilization reagent containing FACSlyse™ + 0.025% Tween 20.

(3) Wash in PBS and stain with MAbs detecting intracellular constituents (e.g. CD69 + IL-2).

Analysis of responding T cells is based on the percentage of T cells which coincidentally express CD69 antigen and cytokines. We have determined that T cell cytokine responses were predominantly restricted to cells expressing the CD69 antigen. A type 1 (TH1) response, for example is typically identified as the percentage of CD3+ T cells co-expressing CD69 and  $\gamma$ -IFN, IL-2, or TNF- $\alpha$ . A type 2 (TH2) response can be characterized as those CD3+ T cells which co-express CD69, and IL-4 or IL-13, for example. The distribution of individual cytokines of type I or type II responses may also be descriptive of a specific kind of cellular response to antigen. For example unique T cell subsets may express TNF- $\alpha$  and  $\gamma$ -IFN in response to different antigens.

Typically the percentage of memory CD4+ T cells co-expressing CD69 and IL-2 or TNF- $\alpha$  in response to a recall antigen such as mumps or CMV antigen is within the range of 0.15% - 3.0% of total CD3+ T lymphocytes. Since this assay procedure measures the frequency of T cell responses to externally processed

antigen - most T cell responses will be class II MHC directed and consequently will be CD4+CD45RO+ T lymphocyte recognitive events. This has been verified by demonstrating that the response to the above mentioned antigens can be blocked by the addition of anti-class II but not anti-class I monoclonal antibodies. Internally processed antigen or class I MHC presented antigen directed responses, however, can be identified by the same set of T cell parameters (i.e. CD69 + cytokine expression).

### BRIEF DESCRIPTION OF THE FIGURES

**Figure 1: Interferon- $\gamma$  response in PBMC and whole blood.** The illustrated flow dot plot histograms compare the  $\gamma$ -interferon response in antigen specific T cells analyzed either from PBMC purified preparations (middle panel) or from whole blood (far right panel) by multiparameter flow cytometry. The first panel represents CMV control antigen response in T cells analyzed from PBMC preparations. Whole blood analysis of the control antigen response was identical (data not shown). These histograms clearly illustrate that detection sensitivity of cytokine expression in T cells analyzed in PBMC and whole blood is similar.

**Figure 2: Mumps antigen response in whole blood.** This figure demonstrates three different T cell cytokine responses to mumps antigen in whole blood. In contrast to figure 1 the T cell response to mumps antigen is primarily characterized by expression of TNF- $\alpha$  with a low frequency of cells expressing  $\gamma$ -IFN. These examples highlight the utility of measuring individual cellular cytokine responses which may vary depending on the source of antigen.

**Figure 3: Drug sensitivity of T cell cytokine expression in response to specific antigen.** This figure serves to show that this assay can be used to examine drug effects on T cell immunity at the single cell level. The upper row of panels indicates the frequency of CD4+ T cells which express various cytokines (IL-1,  $\gamma$ -IFN, and TNF- $\alpha$ ) in response to CMV antigen. The corresponding histograms in the lower panel demonstrate the degree of inhibition of expression by cyclosporine A, a drug used to suppress cellular immune responses in humans. We have also observed both suppressive and enhancing effects of other drugs including hermbimycin, antibodies to accessory molecules and other regulatory cytokines. This figure shows how the described assay can be used to screen drugs which may be targeted to distinct cellular subsets with specific cytokine profiles. This assay may also be used to test the effects of circulating drug in treated patients on specific T cell immunity. The ability to examine cellular responses in whole blood offers the potential to assess drug effects on cellular immunity in a rapid highly informative multiparameter assay.

### **DETAILED DESCRIPTION OF INVENTION**

The instant invention presents a novel method that uses sophisticated, but technically straightforward, flow cytometric method to quantitate and/or phenotypically and/or functionally characterize Ag-specific CD4+ memory/effector human T cells with unparalleled sensitivity in peripheral blood mononuclear cells. Reports in the literature using ELISPOT and limiting dilution type assays suggest precursor frequencies for CMV and mycobacteria reactive T cells in appropriately exposed subjects are on the order of 1/1000 to 1/10,000 mononuclear cell, whereas

the flow cytometric assay described here demonstrates precursor frequencies for these Ags in analogous subject populations to be in the 1/1000 to 2/100 CD4+ T cell range. Although some of these differences may be due to variations in subject populations and Ag preparations, these numbers suggest a significant increase in sensitivity between the flow cytometry-based assay and previous assay types, even after correcting for the different denominators of these assays (e.g. total mononuclear cells for limiting dilution and ELISPOT assays; CD4+ T cells for the flow cytometric assay). We have calculated the sensitivity of the flow cytometric assay for intracellular cytokines using calibrated fluorescent beads and 1:1 molar conjugates of PE to antibody to be on the order of 700 molecules per cell. This value contrast with the picogram quantities of cytokine necessary to detect by ELISA.

The increased sensitivity of the flow cytometric assay is likely due to a combination of factors, including 1) the high sensitivity of fluorescence detection by modern flow cytometers , 2) the highly efficient capture of produced cytokine within the cytoplasm of the secretion-inhibited responding cell (as compared to only partial capture in ELISPOT type assays), 3) the independence of culture (stimulation) conditions and signal detection, allowing these culture conditions to be set up with optimization of response as the only concern, and 4) the relatively short stimulation period, which mitigates against the potential negative effect of activation-induced apoptosis on detection efficiency.

The method described in this invention provides a procedure to detect antigen specific T cell responses in whole blood cultures in under 24 h. It has been found that in whole blood cultures, APCs and T cells form adhesive complexes as a consequence of specific antigen activation which prevented detection of responding

T cells in cell suspension assays like flow cytometry. The essential modification described herein defines methods for quantitatively recovering small numbers of T cells(<0.1%) from whole blood cultures for analysis from the adhesive antigen-presenting cell (APC) cellular complexes.

5

Increased sensitivity is not the only advantage of the flow cytometry based assay of this invention. In addition to simplicity, which makes it amenable to routine clinical laboratory testing, this assay also has the significant advantages of being able to assess 1) Ag-specific responses in phenotypically defined T cell subsets, and 2) the synthesis of multiple cytokines in a single Ag-responsive T cell. With respect to the subset delineation, we have used the multiparameter capabilities of this assay to precisely define Ag-specific responses within T cell subsets defined by CD4 vs. CD8, TCR classes, homing receptors, and markers such as CD45RA or CD27. Such subset assignment of the Ag-specific responses allows the separate evaluation of Ag-specific responses in distinct T cell subsets that are potentially independently regulated (e.g. CD4+ vs. CD8+ T cells), and greatly increases the precision and reproducibility of these assays by automatically correcting for potentially confounding changes in the overall representation of a particular T cell subset within the overall PBMC population. Obviously, this capability is particularly germane for the study of CD4+ T cell effector frequencies in HIV+ subjects, whose CD4 counts may vary tremendously.

With regard to the ability to assess multiple cytokines per cell, this assay is the only method described to date that allows the delineation of cytokine secretion patterns of freshly-isolated T cells on an Ag-specific, single cell basis. Although the concept of cytokine synthesis heterogeneity among memory/effector T cell populations (e.g. the TH1 vs. TH2 paradigm) has been well established, it is unclear



the extent to which the described cytokine synthesis phenotypes of cloned T cells represent the true spectrum of memory/effector T cell cytokine synthesis function *in vitro*. Cloned T cells, by definition, spend weeks to months *in vitro* prior to analysis -- during which time functional biases may be introduced by clonal selection and/or inadvertent regulation of effector cell differentiation. The feasibility of determining cytokine synthesis patterns of Ag-specific T cells taken directly from the blood, and stimulated *in vitro* for only 6 hours. While the predominant pattern of cytokine production of CMV-specific CD4<sup>+</sup> T cells from normal subjects is consistent with a classical TH1 phenotype ( $\gamma$ -IFN<sup>+</sup>/IL-4<sup>-</sup>), variably-sized subsets of CD4<sup>+</sup> T cells synthesizing all or most of the possible combinations of these cytokines on a single cell basis have also been observed.

## EXAMPLES

The following examples illustrate certain preferred embodiments of the invention but are not intended to be illustrative of all embodiments.

### Example 1

#### **Cell preparation and antigen stimulation:**

Heparinized whole blood samples were collected from normal or diseased donors using heparin Vacutainer<sup>TM</sup> blood collection tubes (Becton Dickinson VACUTAINER Systems, Franklin Lakes, NJ). One ml aliquots of whole blood were dispensed into 16X125 polypropylene tissue culture tubes (Corning Costar Corporation, Cambridge, MA).

Appropriately titred specific or control Ag preparations (60µl/ml) and in most instances (except as noted), 4 µg anti-CD28 MAbs were added to 1ml whole blood aliquots. Cultures were incubated at a 5° slant at 37°C in a humidified 5% CO<sub>2</sub> atmosphere for 1 hour with slight agitation to improve APC/T cell interaction, and an additional final 5 hours including a final concentration of 10µg/ml of Brefeldin A (a relatively non-toxic, but potent, inhibitor of intracellular transport which prevents secretion of any produced cytokines. 20mM Na<sub>2</sub>EDTA for a final concentration of 2mM was then added directly to the whole blood cultures. Samples were vortexed and incubated for 15 min. @ RT, washed in PBS with vigorous pipetting. Whole blood samples were subsequently RBC lysed and fixed with FACSlysing™ solution (Becton Dickinson) and washed in PBS. plus 1% BSA prior to resuspension @ 5 x 10<sup>6</sup> cells/ml in freezing medium (10% dimethyl sulfoxide in dPBS with 1% BSA). Finally, the cells were frozen in 2ml Wheaton vials at -70°C in a Nalgene (Rochester, NY) freezing chamber.

## **Example 2**

### **Immunofluorescent staining:**

Cell preparations frozen as described above were rapidly thawed in a 37°C water bath and then washed once with cold dPBS prior to resuspension in fixation/permeabilization solution (Becton Dickinson Immunocytometry Systems, San Jose, CA) at 1 x 10<sup>6</sup> cells/500ul, and incubation for 10 minutes at room temperature in the dark. These cells were washed twice with dPBS with BSA and sodium azide, and then were stained protected from light with directly conjugated MAbs for 30 minutes.. In some experiments, the freezing step was omitted, and cells freshly harvested after Ag activation were cell surface stained first, and then fixed/permeabilized/washed, and then stained for intracytoplasmic Ags. After

staining, the cells were washed, refixed in 1% paraformaldehyde in dPBS, and then kept protected from light at 4° C until analysis on the flow cytometer.

### **Example 3**

#### **5 Flow Cytometric Analysis:**

Up to six parameter analysis was performed on a modified FACSort™ (Becton Dickinson) flow cytometer equipped with a second 532 nm line diode laser (BDIS) using FITC, phycoerythrin(PE), and peridinin chlorophyll protein (PerCP), and allophycocyanin (APC) as the 4 fluorescent parameters, using methods of cytometer set up and data acquisition which are well known in the art. Otherwise 5 parameter analysis using up to three fluorescent parameters was conducted on a standard FACScan™ flow cytometer. For a typical analysis, 50,000 events were acquired, gated on CD4 expression and a light scatter gate designed to include only viable lymphocytes (most files required "fine-tuning" of gating during analysis on these same parameters, leaving 48,000 events for the final profiles). In some analyses, additional "live" gating based on a fluorescent parameter (for example, CD69 reactivity) was performed to enhance the sampling of small populations. List mode multiparameter data files (each file with forward scatter, orthogonal scatter, and 3-4 fluorescent parameters) were analyzed using the PA1NT~A~GATEPlus™ program (BDIS). Isotype-matched negative control reagents were used to verify the staining specificity of experimental antibodies, and as a guide for setting markers to delineate "positive" and "negative" populations.

#### **Example 4**

#### **Establishment of a Flow Cytometric Assay of Ag-Specific T Cell Cytokine Responses:**

5           Previous methods for detecting and quantitating T cell cytokine responses to polyclonal mitogens (e.g. phorbol ester plus ionomycin) and superAgs (the staphylococcal enterotoxin superAgs SEA and SEB) using a protocol based on intracytoplasmic staining of cytokine in short-term activated, secretion-inhibited T cells and multiparameter flow cytometric analysis has been previously known. It  
10 was initially anticipated that, with slight modifications, this approach would identify and quantitate T cells producing cytokine in response to nominal Ags as well. However, Ag responses differed from mitogen and superAg responses in several key areas. First, the geometry of the T cell/accessory cell interaction was critical for Ag responses; maximal responses were observed in slant tubes that allowed close  
15 proximity of T cells and accessory cells, but still allowed adequate media access to responding cells. Second, responses were maximized when Brefeldin A was omitted from the initial hour of interaction (likely to allow optimal antigen processing), and when exogenous costimulation was provided. Third, precise detection of responding T cells was enhanced with inclusion of CD69 assessment in  
20 the multiparameter protocol. CD69 is upregulated on activated T cells prior to cytokine production, and thus allows more definitive "clustering" of the true responding fraction. Background staining, when present, is often only present in the CD69 negative fraction, and thus can be excluded from consideration. CD69 can also be used as a live gating parameter to enhance collection of the cytokine  
25 producing cells, yet still allow back calculation of responding cell frequency to overall CD4+ population. Finally, because of the relatively small size of the Ag-

specific populations, accurate assessment of these responses required the routine collection and analysis of at least 50,000 events per determination.

In order to allow antigen specific reactions to occur in whole blood environments additional modifications were required. Key elements in these modifications were: (1) Slight agitation during the first hour of incubation of whole blood cultures with antigen to improve APC interactions with specific T cells, (2) use of polypropylene culture tubes to more easily displace T cells from adherent APC complexes, (3) direct addition of EDTA to whole blood cultures to optimize recovery of antigen-activated T cells. Figure 1 compares frequencies of CD4<sup>+</sup> T cells responding to CMV antigen. This figure illustrates similar T cell response behavior in both PBMC and whole blood environments (e.g. when samples are obtained from normal donors). Figure 2 demonstrates that T cell responses to other cognate antigens such as mumps virus can also be measured using whole blood culture methodology. In contrast to the CMV response example 1, T cell responses to mumps in this individual exhibits a proportionally higher frequency of T cells expressing TNF- $\alpha$  and a smaller frequency of CD4<sup>+</sup> T cells expressing  $\gamma$ -IFN.

## **Example 5**

### **Assessment of drug effects on antigen specific reactions in whole blood:**

Concentrations (typically less than 1 $\mu$ g of immunosuppressive drug (e.g. Cyclosporine A, Herbimycin A) were added to initial whole blood cultures activated by antigens as described in example 1. Frequencies of CD4<sup>+</sup> T cells responding to various antigen in the presence or absence of drug were compared as determined by flow cytometry as detailed in example 2 and example 3. Figure 3 (example 3)

shows the suppressive effects of the addition of 1µg of cyclosporine A to antigen specific T cell responses in whole blood cultures. In this example, 75% inhibition of IL-2, 87% of γ-IFN and 92% of TNF-α T cell responses to CMV were blocked by cyclosporine A. This example shows how the whole blood assay can be used to assess immunosuppressive or augmenting effects of immunomodulatory drugs. A corollary to this assay -is that effects of in vivo concentrations of drugs may be also be assessed. Since the assay is performed in the presence of autologous plasma, T cell responses can be measured in the presence of pharmacological concentrations of administered drugs.

It is apparent that many modifications and variations of this invention as hereinabove set forth may be made without departing from the scope and the spirit thereof. The specific embodiments are given by way of example only, and the invention is limited only by the terms of the appended claims.

### **WHAT IS CLAIMED IS:**

1. A method for determining the antigen specific activation of T cells from a human or animal subject comprising:
  - (i) culturing a sample of whole blood of said human or animal with an antigen specific stimulus; and
  - (ii) analyzing an aliquot of said cultured sample for expression of one or more intracellular cytokines and/or one or more early activation antigens in combination with one or more T cell markers, by means of multi-parameter flow cytometry.
2. The method of Claim 1 further comprising lysing the red blood cells and washing to remove the debris prior to the measuring step.
3. The method of Claim 1 further comprising providing exogenous costimulation of said cells by adding a costimulus to said sample contemporaneous with the addition of said antigen specific stimulus.
4. The method of Claim 3 wherein the costimulus is selected from the group consisting of ligands or antibodies to CD28, CD40, CD86, and CD118.
5. The method of Claim 1 wherein the antigen specific stimulus is selected from the group consisting of alloantigens, viral antigens, autoantigens, peptides, and proteins.

6. The method of Claim 5 wherein the antigen specific stimulus is selected from the group consisting of cytomegalovirus antigen, measles virus antigen, Mycobacterium tuberculosis antigen, and mumps virus antigen.
7. The method of Claim 1 wherein said intracellular cytokines are selected from the group consisting of  $\gamma$ -IFN, IL-2, IL-4, and TNF- $\alpha$ .
8. The method of Claim 1 wherein early activation antigen is CD69.
9. The method of Claim 1 wherein the sample is cultured with the antigen specific stimulus for a period of 24 hours or less.
10. The method of Claim 1 wherein the sample is cultured with the antigen specific stimulus for a period of 6 hours or less.
11. The method of Claim 1 further comprising culturing the sample with an immunosuppressive or immuno augmenting drug, simultaneous with the culturing with the antigen specific stimulus.
12. The method of Claim 1 wherein the sample of whole blood is obtained from a human or animal treated with an immunosuppressive or an immunoaugmenting drug.
13. The method of Claim 1 further comprising adding a reagent which blocks secretion of intracellular cytokines during the culturing of the sample with the antigen specific stimulus.
14. The method of Claim 13 wherein the reagent is Brefeldin A.



15. The method of Claim 13 wherein the reagent is added after the culturing has proceeded for at least 1 hour.

16. The method of Claim 1 further comprising adding a cationic chelating agent after the culturing is complete but prior to the measuring.


17. The method of Claim 16 wherein the cationic chelating agent is EDTA or EGTA.

18. The method of Claim 1 wherein the lysing of the red blood cells is accomplished using a lysing agent.

## **ABSTRACT**

This invention comprises a novel approach to the assessment of antigen-specific T cells that quantitates and characterizes these cells with unprecedented clarity, and importantly, because it is performed in whole blood, is amenable to routine use in the clinical immunology laboratory. The methodology offers an improved flow cytometric intracellular cytokine assay in whole blood that can simultaneously measure multiple T cell subsets expressing multiple cytokines from a single whole blood culture. Evaluation of whole blood antigen specific cytokine responses has the important advantage of assessing T cell activation in the presence of ALL types of MHC autologous antigen presenting cells present in the native sample. It also has the advantage of enabling a culture system (whole blood) which can reflect effects of systemic environments (i.e. drug augmentation or suppression) on T cell responses to specific stimuli including antigen. by either culturing in the presence of such drug or analyzing the blood of a human or animal receiving such drug.


## IFN $\gamma$ Response in PBMC and Whole Blood (CD4 gated)

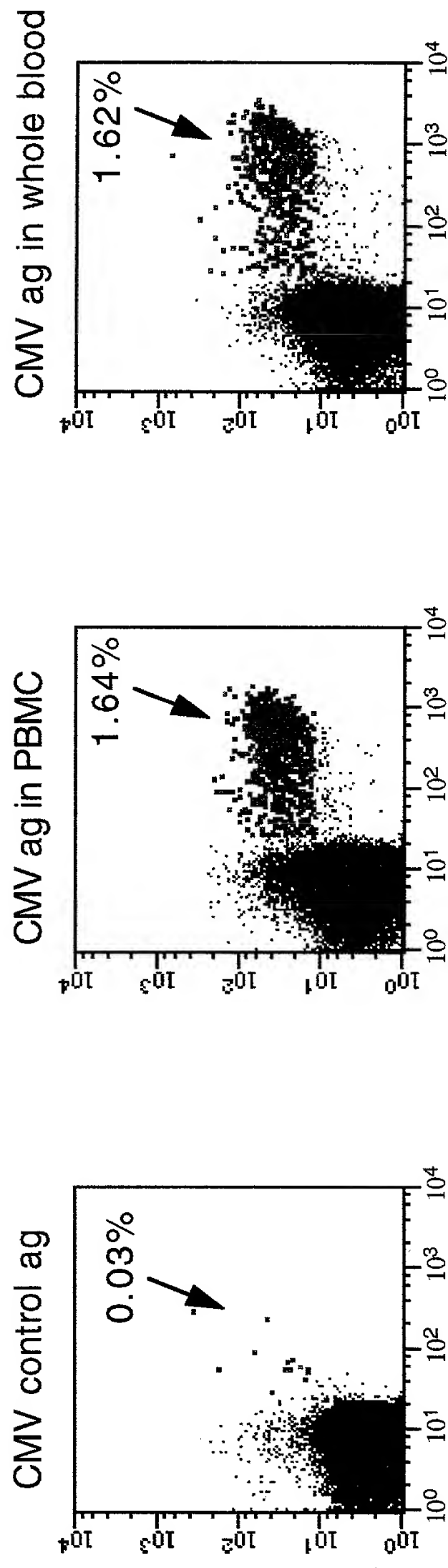


$\alpha$ IFN $\gamma$  FITC vs. CD69 PE

FIGURE 1


## Mumps Antigen Response in Whole Blood (CD4 gated)

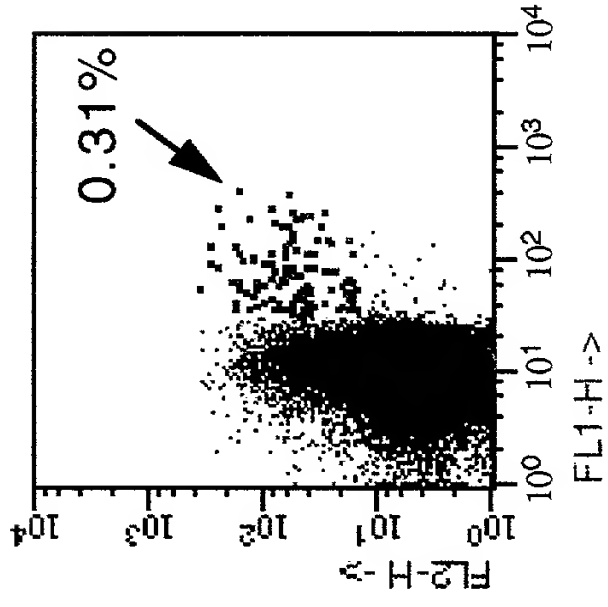

αIL-2 FITC/CD69 PE

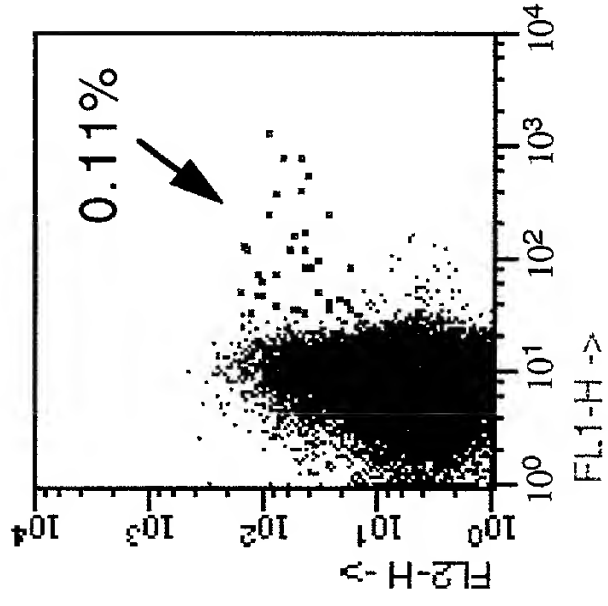

αIFN $\gamma$  FITC/CD69 PE

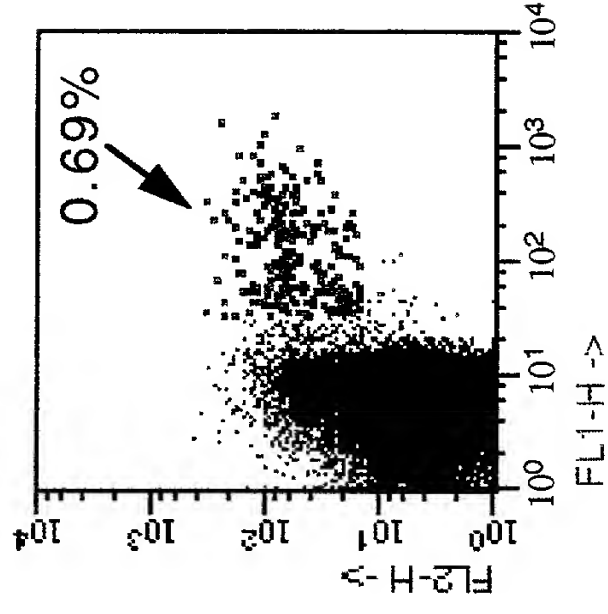

αTNF $\alpha$  FITC/CD69 PE

FIGURE 2

# Cyclosporin A Inhibition of CMV ag Response (CD4 gated)

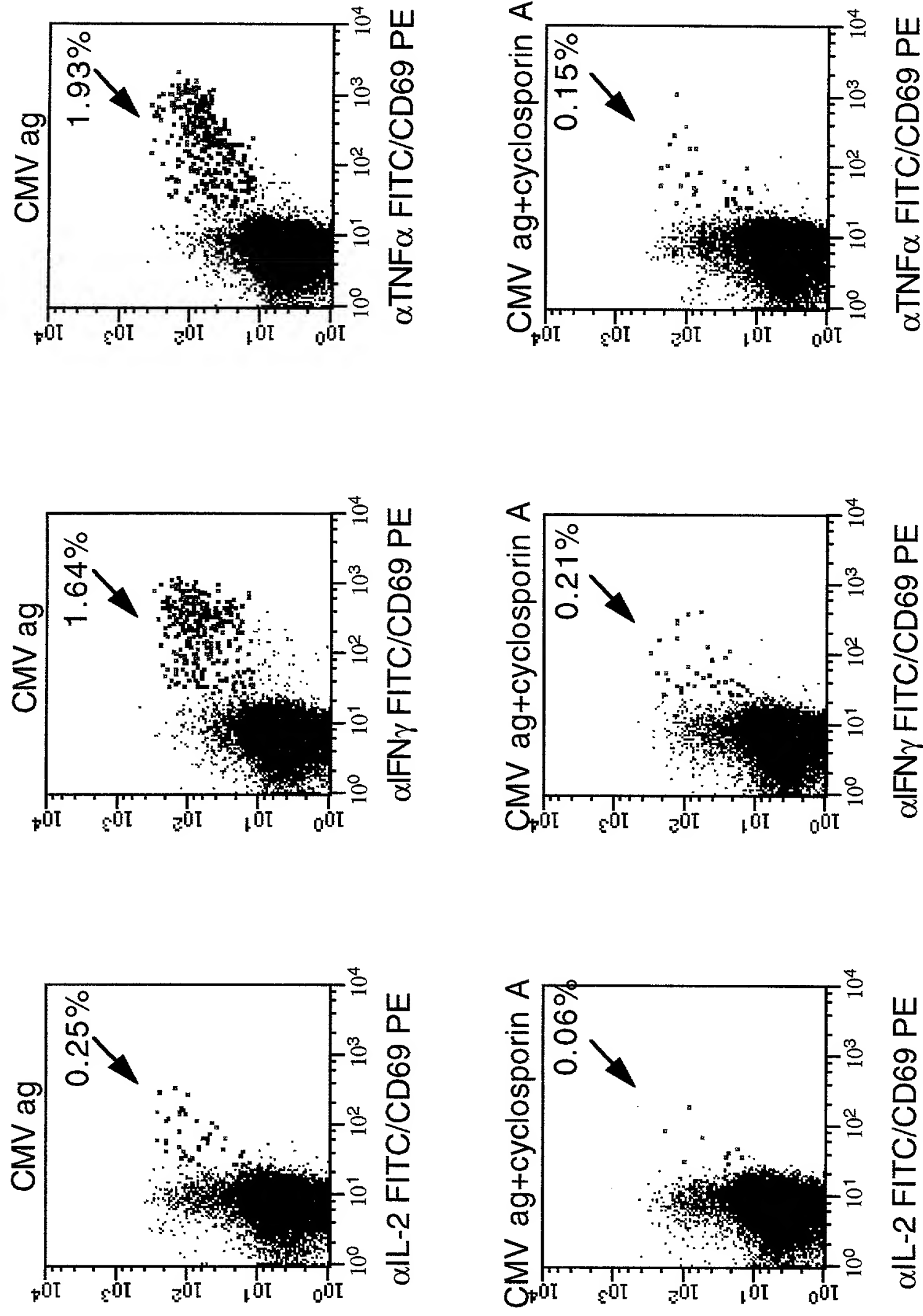


FIGURE 3